

# Capp *Rondo*


TUBE ROLLER

**Capp *Rondo* Tube Roller** ensures efficient mixing of blood samples, viscous substances and liquid-solid suspensions. The programmable pulse mode feature with Run-Pause-Run motion and variable time setting (30 to 99 seconds) give the user additional flexibility and improve the mixing process. In addition to the rolling motion, the rollers have a rocking movement which ensures optimal mixing throughout. Equipped with a brushless DC motor, **Capp *Rondo* Tube Roller** provides quiet and maintenance-free operation, highly suitable for extended runs.

## Highlights:

### *Improved mixing performance:*

- In addition to the rolling motion, the rollers of the **Capp *Rondo* Tube Roller** rock by lowering and raising, which improves the mixing process.
- Programmable momentary pulse mode ensures throughout optimal mixing.

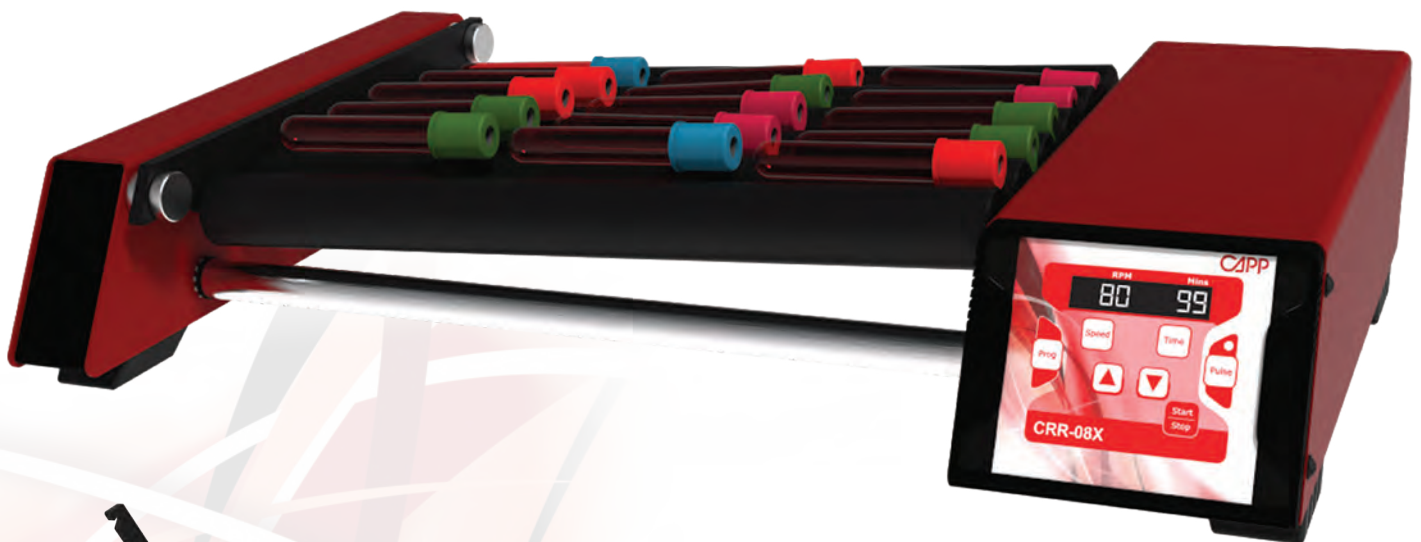

Removable Rollers



Stackable System



**Capp Rondo Tube Roller** can accommodate up to 15 blood tubes of 7mL. The user interface panel is highly intuitive and displays both the time and the speed.

**Capp Rondo Tube Roller** operates at the speed range from 10 up to 80 rpm. The user can program up to 9 customized protocols.

Focusing on increased user comfort and the low contamination risk, **Capp Rondo Tube Roller** operates with easily removable and autoclavable rollers. The stackable design and relatively small footprint will optimize your bench space. **Capp Rondo Tube Roller** is also suitable for use in incubators with temperatures up to 40°C and humidity up to 80%.



#### Highlights:

- Suitable for incubators and highly humid environments.
- Tool free detachable and autoclavable rollers for increased user comfort and lower contamination risk.
- By adjusting the rollers space, you can accommodate multiple sizes of tubes and bottles.
- Stackable design with a small footprint for optimal bench space use.

Specifications:	Capp Rondo Tube Roller
Motor type	Brushless DC motor
Motion	Rocking and rolling
Number of tubes	up to 15 x 7ml blood tubes
Speed range	10 to 80 rpm
Time range	1 to 99 mins
Programmable range	User can set and save up to 9 user-defined programs (protocols)
Operation type	Continuous and timed operation
Removable rollers	Yes, tool free removal and fitment
Pulse mode	Yes (Run-Pause-Run mode)
Tilt angle	+/-3°
Protection class	IP21
Weight	6kg
Dimensions (W x D x H)	502mm x 92.5mm x 301mm
Warranty	2 years

## Ordering Info

Cat No.	Description
CRR-08X	CappRondo Tube Roller 80 rpm